



Polyurethane Enteral Feeding Tube

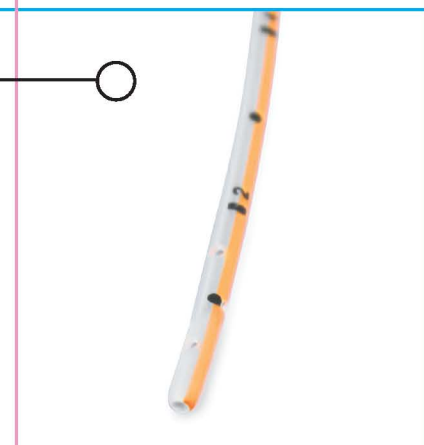

The open end design, with precision drilled side holes, enhances fluid flow and decreases potential feeding stagnation and mucosal irritation.

The hub is clearly marked with "ORAL" and French size to readily identify the tube as an enteral feeding tube.

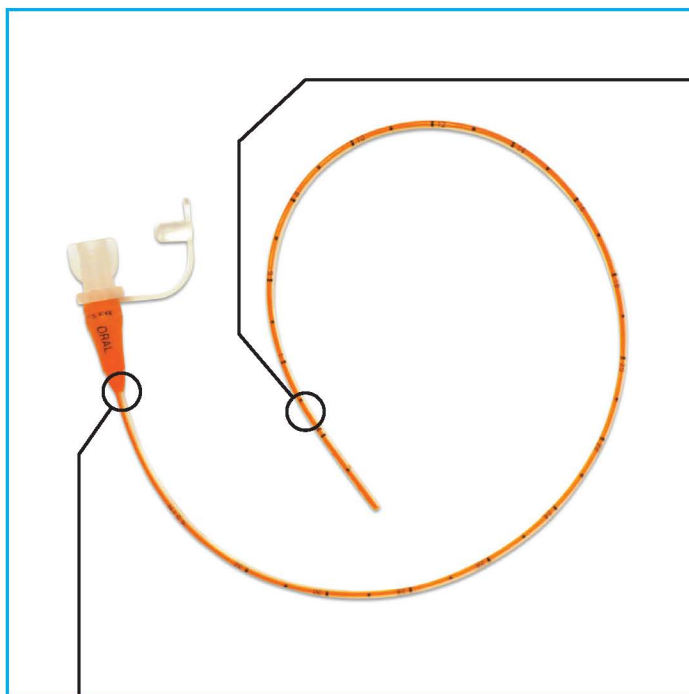

Polyurethane Enteral Feeding Tube

- Unique polyurethane tube with silicone overlaid hub designed to provide *secure* enteral only connections
- Labeling supports ISMP recommendations that all distal ports and tubing of access lines clearly indicate what is the use of the port/line
- Bold, indelible markings every centimeter
- Color coded for patient safety
- Polyurethane's tensile strength allows thin wall construction designed for increased flow rates and decreased clogging with proper maintenance
- Latex and DEHP free
- Integrated closure plug for easy cleaning
- Radiopaque stripe for confirmation of tip placement by x-ray



Footprint
medical inc.

210.226.2600

www.footprintmed.com

Ordering Information for:

Polyurethane Feeding Tubes



Polyurethane 40cm Enteral Feeding Tubes

P1FTS4.0-EO	Polyurethane Feeding Tubes	4.0 Fr.	10/cs
P1FTS5.0-EO	Polyurethane Feeding Tubes	5.0 Fr.	10/cs
P1FTS6.5-EO	Polyurethane Feeding Tubes	6.5 Fr.	10/cs
P1FTS8.0-EO	Polyurethane Feeding Tubes	8.0 Fr.	10/cs

Polyurethane 60cm Enteral Feeding Tubes

P1FTM4.0-EO	Polyurethane Feeding Tubes	4.0 Fr.	10/cs
P1FTM5.0-EO	Polyurethane Feeding Tubes	5.0 Fr.	10/cs
P1FTM6.5-EO	Polyurethane Feeding Tubes	6.5 Fr.	10/cs
P1FTM8.0-EO	Polyurethane Feeding Tubes	8.0 Fr.	10/cs

Polyurethane 90cm Enteral Feeding Tubes

P1FTL5.0-EO	Polyurethane Feeding Tubes	5.0 Fr.	10/cs
P1FTL6.5-EO	Polyurethane Feeding Tubes	6.5 Fr.	10/cs
P1FTL8.0-EO	Polyurethane Feeding Tubes	8.0 Fr.	10/cs
P1FTL10.0-EO	Polyurethane Feeding Tubes (110cm)	10.0 Fr.	10/cs

"Helping tiny babies make their own footprints in the world..."

Distributed By:



Phone: 847.891.6210
Fax: 847.891.6215
www.accessmedicalsupply.com

Access Medical Supply, Inc.



12727 Cimarron Path
San Antonio, TX 78249
210.226.2600
www.footprintmed.com